



Surgical Light

AIM-200 OR

Versatile, flexible and high performing.

PHILIPS burton

#### Surgical Light

# AIM-200 OR


The right light for surgical procedures.

Burton's AIM-200 OR provides versatility, amazing illumination, and ease of use. The rotating center handle makes adjusting the light size and pattern to fit the working area quick and convenient. This light's superior arm mechanism is easy to position and will not drift. The AIM-200 OR is available with a floor, wall, single-, or double-ceiling mount which makes it a practical choice for hospitals, ambulatory surgery centers and veterinarian offices.

#### Features

- 90,000 lux (8361 fc) at 1 meter
- 3700 K color temperature
- CRI (Color Rendering Index) of 96
- Larger and adjustable light pattern
- Unlimited rotation and maximum reach of 63" (160 cm)
- Light field diameter: 10" 12" (25 cm 30 cm)
- Depth of illumination: 47" (119 cm)
- Adjustable focusing by rotating center handle
- UL60601-1 / IEC 60601-1 / IEC 60601-1-2 / IEC 60601-2-41; CAN/CSA C22.2 601-1 M90 certified
- CE marked
- 5-year limited warranty
- Assembled in USA



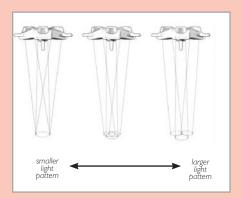


#### Larger and Adjustable Light Pattern

The multiple reflectors of the AIM-200 OR can be tilted by rotating the sterilizable light handle to adjust the size and pattern of the light to the working area.



Light Pattern Free of Cast Shadows

The AIM-200 OR has a unique Y-Shaped design, which facilitates placement of the lighthead around the practitioner's head minimizing the obstruction of light into the surgical field.

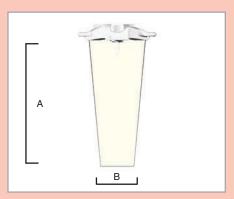

#### Depth of Field / Beam Diameter

		cm i	n
А	Depth of Field	119	47
В	Beam Diameter D10	25	10




#### Ordering Information

120V Model	230/240V Model	230V SP Model*	100V Model	Description
A200FL	A203FL	A202FL25	A201FL	Floorstand
A200SC	A203SC	A202SC25	A201SC	Single Ceiling Mount
A200DC	A203DC	A202DC25	A201DC	Double Ceiling Mount
A200W	A203W	A202W25	A201W	Wall Mount
Accessories				
1017111PK				AIM-200 OR Replacement Bulbs, 3/Box
1017040				Sterilizable / Autoclavable Handle
0008100PK				Light Handle Covers, Sterile, 25/Box
1017084				42" Drop Tube

\*SP Model refers to Spanish Language Manuals.

## Unlimited Rotation and Maximum Reach

The AIM-200 OR's 360° limitless arm-and-mounting-system-rotation around vertical axes on single ceiling models is a key feature previously available only on more expensive OR lights. The arm system also provides a maximum reach of 63" (160 cm) and diameter of 10' (305 cm), allowing it to illuminate the patient from head to toe with no drift.

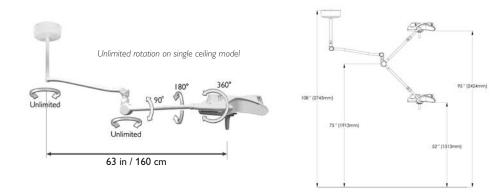

### Summary of technical data

90,000 lux (8361 fc) at 1 meter		
3700 K		
96		
20" (51 cm)		
10" - 12" (25 cm - 30 cm)		
47" (119 cm)		
Adjustable by rotating center handle		
50W IRC, 12V, bi-pin		
3		
1,400 hours		
60" (152 cm) max.		
41" (104 cm) vertical movement		
8 ft./12 ft. (244 cm/366 cm)		
105 Watts		
Floor version 35 lbs. (15.9 kg), single ceiling version 34 lbs. (15.4 kg)		
Yes		
UL60601-1 / IEC 60601-1 / IEC 60601-1-2 / IEC 60601-2-41		
CAN/CSA C22.2 601-1 M90		
CE marked		
5-year parts and labor		
Assembled in USA		

Note: The optical data are at nominal values based on measurements done according to the UL/IEC performance standards at 120V \* With extension arm and spring arm extended

© 2016 Koninklijke Philips N.V. All rights reserved. Philips reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Phone (800) 444-9909, Fax (800) 765-1770



Philips Burton 11500 Melrose Avenue Franklin Park, IL 60131 Burtonmedical.com

